


New Walser[®] Tooth Matrix! 4-surface fillings in half the time

At the world's leading trade fair in Cologne, this innovation of **hundreds of thousands of products was elected one of the highlights of the IDS 2019** and immediately ordered by dentists worldwide.

At the request of the dentists who already work with Walser[®] matrices in more than 120 countries, the innovative dental matrix has been developed from the practice for the practice.

This Walser[®] matrix with automatic and adjustable adaptation can be adapted e.g. palatally depending on the extent of the tooth shape to be built up.

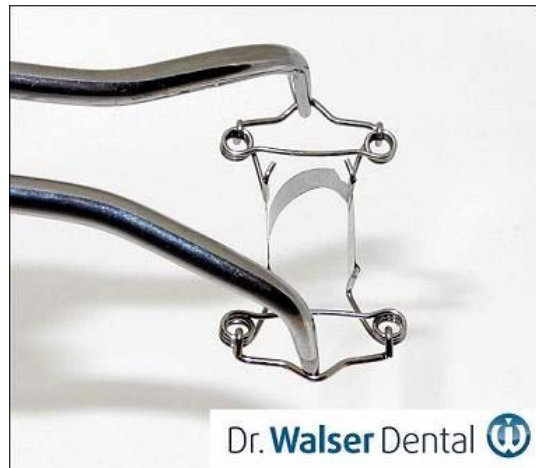

Order now

The application is simple, the time required for the reconstruction is almost halved:

1. While tensioning the matrix, at the same time press down the raised band with your finger so that it is the same height as the contacting band.
2. By stretching the matrix above the molar, you control the parallel displacement of the bands yourself. You can rotate the matrix on both the palatal and buccal side, the same applies to the lower jaw.
3. When releasing the matrix on the tooth, the band automatically slides in the interdental direction. You can use an instrument (such as a plugger) to control yourself how far the band is pushed into the interdental space.



Tens of thousands of dentists in more than 120 countries worldwide work with WALSER[®] Matrices and that number increases every day...

Dr. Walser Dental GmbH

E-Mail: info@walser-dental.com

Geschäftsführer/Managing Director:

Fritz-Reichle-Ring 18

Phone +49 77 32 – 33 00

Senator h.c. Gerhard R. Daiger

78315 RADOLFZELL/GERMANY

Fax +49 77 32 – 5 72 23

VAT-ID: DE142779045

www.walser-dental.com

HRB: 550113, Freiburg i. Br.